

MAGIO MS-601F Refrigerated / heating circulator

As with all circulators from the MAGIO range, the refrigerated circulators stand out thanks to their premium quality, high performance and intuitive operation. The devices offer extra strong pressure and suction pumps, thus fulfilling the highest demands for temperature control of external applications. Whether in basic research, material testing or technical systems – the MAGIO refrigerated circulators offer high-tech solutions for high customer requirements.

High resolution TFT touch display

The modern TFT touch display gives you all the important information at a glance. Three large, predefined main screens clearly display data and graphics with various application priorities. Menu navigation is self-explanatory, arranged by relevance to daily operations and easy to operate with the touch of a finger. The in-built help function provides detailed support in case of additional questions.






Product features

- · Ideal for demanding external applications
- Simple control of complex applications
- Continuously adjustable, extremely powerful pressure / suction pump
- Flow rate 16 ... 31 I / min, pressure 0.24 ... 0.92 bar, suction 0.03 ... 0.4 bar
- Large, high-resolution TFT touch display with multilingual user interface
- · Parts being in contact with the medium made of stainless steel
- Integrated programmer
- Integrated external Pt100 connection
- USB connection
- RS232 interface for online communication
- Ethernet
- analog interfaces (accessory)
- Class III (FL) according to DIN 12876-1
- Connections for solenoid valve
- Integrated pump connection M16×1



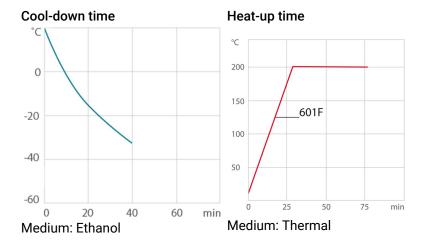

Performance values

230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)		
Heating capacity kW	2	
Viscosity max. cSt	70	
Pump capacity flow rate I/min	16 31	
Pump capacity flow pressure psi	3.5 13.3	
Maximum suction psi	-0.45.8	
Power consumption A	16	

0000705.00						
Order No.		9032705.33				
Cooling capacity	(Ethanol)					
°C	20	10	0	-10	-20	-30
kW ¹	0.6	0.52	0.44	0.27	0.16	0.04
Refrigerant stag	e 1					
Refrigerant		R449A				
Filling weight g		150				
Global Warming R449A	Potential for	1397				
Carbon dioxide	equivalent t	0.20955				

¹ Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

Technical data

Available voltage ve	ersions	Cooling	
Order No.	9 032 705	Cooling of compressor	1-stage Air
Available voltage versions:			
9032705.01	100V/50-60Hz (Nema N5-15 Plug) (R452A)		
9032705.02	115V/60Hz (Nema N5-15 Plug) (R449A)		
9032705.05	200-230V/50-60Hz (CH Plug Type SEV 1011) (R449A)		



9032705.04	200-230V/50-60Hz (UK Plug Type BS1363A) (R449A)
9032705.33	200-230V/50-60Hz (Schuko Plug - CEE 7/4 Plug Type F) (R449A)
9032705.33.chn	200-230V/50-60Hz (CN Plug) (R449A)

Bath		Other	
Bath tank	Stainless steel	Classification	Classification III (FL)
Bath cover	integrated	IP Code	IP 20
Usable bath opening in. (W x L / D)	8.7 x 5.9 / 7.9	Pump function	Pressure Suction Pump
		Pump type	Immersion Pump
		User Interface Language	Chinese, Czech, Dutch, English, French, German, Italian, Japanese, Korean, Polish, Portuguese, Russian, Spanish, Turkish

Ethernet, Modbus, RS232, RS485, Stakei, USB
integrated
8x60 steps
ICC
10 Point Calibration
7" TFT Touchscreen
Touchscreen
00:00 99:59

Dimensions and volumes	
Weight lbs	91.5
Dimensions in. $(W \times L \times H)$	13 x 18.5 x 29.1
Filling volume l	8 10
Pump connections	M16x1 male

Temperature values	
Setting the resolution of the temperature display °C	0.01
Working temperature range °C	-35 +200
Temperature stability °C	±0.01
Ambient temperature °C	+5 +40
Setting the resolution of the temperature display °C	0.01

Included in delivery 2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male)

All Benefits



Intelligent temperature control.

Intelligent cascade control - automatic and self-optimizing adaptation of the PID control parameters with external stability of +/- $0.05\,^{\circ}$ C.



Many interfaces.

Straight-forward remote control, data management, and integration into process structures. USB, Ethernet, RS232, SD card, and alarm off are permanently integrated. Further interfaces available as accessories.





Touch display. Perfect operation.

With the touch display, the user always has an overview of all values and functions. The intuitive and multilingual menu structure enables perfect control.



Maximum safety.

Classification III according to DIN12876-1 enables safe operation, even with flammable fluids.

Automatic switch-off in the event of high temperature or low liquid level.



Space saving. Free up space.

Place your JULABO Circulator right next to an application, another unit, or wall. That saves space. This is made possible by eliminating vents and connections on the sides.



Multi-lingual.

Operation in multiple languages.



Programmer. Integrated.

The integrated internal programmer makes it possible to automatically run temperature time profiles.



Analog I/O.

Analog interfaces for integration into process control systems (optional).



Temperature. Under control.

External Pt100 sensor connection for precise measurement and control directly in the external application.



Fill level. Monitored.

Fill level indicator on the display for heat-transfer liquid.



Process stability.

Early warning - visual and acoustic - of critical states increases process stability.



Process. Under control.

Full control of the dynamic, access to all important control parameters for individual process optimization.



Stable. Mobile.



Energy-saving.

The high-quality insulation of all relevant components saves energy.



Everything made of stainless steel.

Quality and material compatibility at the highest level. All parts in contact with the medium are entirely made of stainless steel.



Most powerful pump.

The integrated pressure/suction pump with performance values of 0.9 bar and -0.4 bar is the most powerful in its class and continuously adjustable.



Wide range.

Refrigerated and heating circulator in various combinations, circulator in various sizes.

Maximum flexibility through a large selection of accessories.



Condensation protection.

Superb design solution. Integrated ventilation directs air over the bath lid and minimizes condensation.



Connection. Easy.

Inclined pump connections (M16×1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



Highest measuring accuracy

'Absolute Temperature Calibration' for manual compensation of a temperature difference, 10-point calibration



100% Checked.

100% testing. 100% quality. Each JULABO Circulator undergoes thorough quality testing before leaving the factory.



Green technology.

Development consistently applied environmentally friendly materials and technologies.





JULABO. Quality.

Highest standards of quality for a long product life.



Quick start.

Individual JULABO consultation and comprehensive manuals at your disposal.



Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.com.



Quick support

If an error occurs, the integrated Black-Box function permits fast diagnosis by the JULABO service team